



A Dräger and Siemens Company

Field Service Procedure

Part Number: SP00284

Rev: -

Date: 15 September 2003

© 2003 Draeger Medical, Inc.

AGS Scavenger Mounting

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

AGS Scavenger Mounting

Purpose: The purpose of this document is to provide instruction for the installation of the AGS Scavenger Mount on the Fabius GS and Tiro, mounting of the Scavenger, and installation of the Draeger Clamp.


Note: The Mounting Bracket is pre-assembled for left side attachment on the Fabius GS and must be reconfigured for other mounting options. Refer to the photos below for all Mounting Bracket configurations and follow the applicable steps. The Draeger Clamp is for use with the right mount COSY on Fabius Tiro.

FABIUS GS ASSEMBLY INSTRUCTIONS :

Configuring the Mounting Bracket for the Application

1. Using a Philips screwdriver, remove the Slide (item 1) by removing 4 screws from the existing configuration and reattach the Slide in the new configuration as shown below.
2. Using the 5/32" hex key (provided), remove the two (2) mounting posts (item 2) from the existing configuration and reattach them in the new configuration as shown below. Ensure that the mounting posts are located in the upper two mounting holes on the bracket.



Mounting Bracket Configurations for Fabius GS

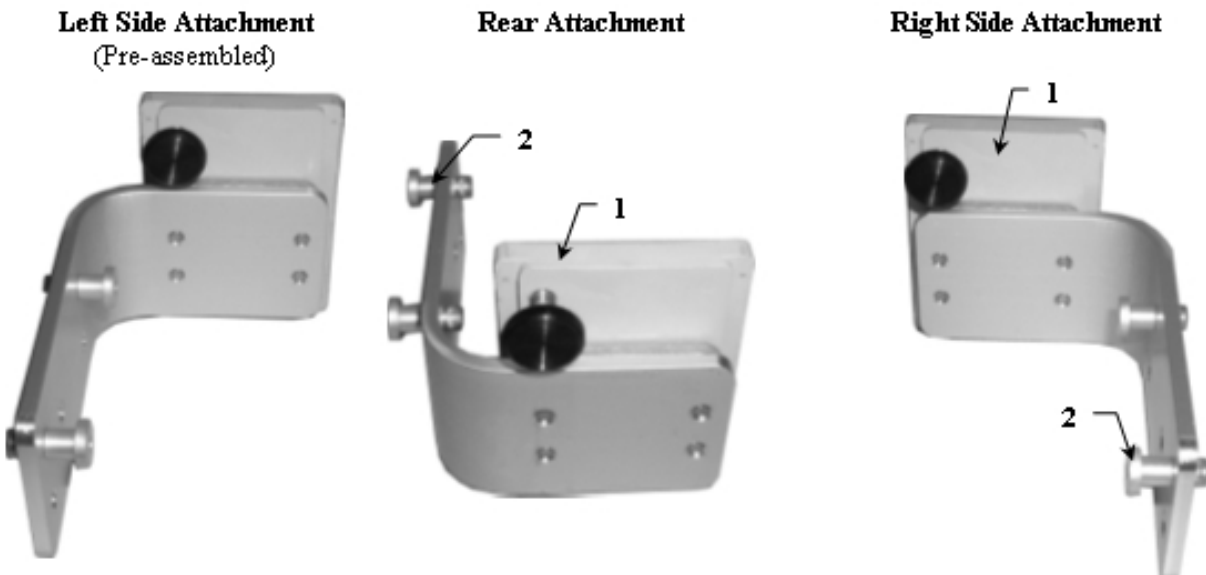

SP00284
9/18/03

Attaching the Scavenger to the Mounting Bracket for Fabius GS

1. Guide the two (2) mounting posts into the mounting slots (item 1) on the rear of the Scavenger.
2. Using the 5/32" hex key (provided), secure Scavenger to Mounting Bracket with the M5 x 18mm socket head cap screw (SHCS) (item 2).

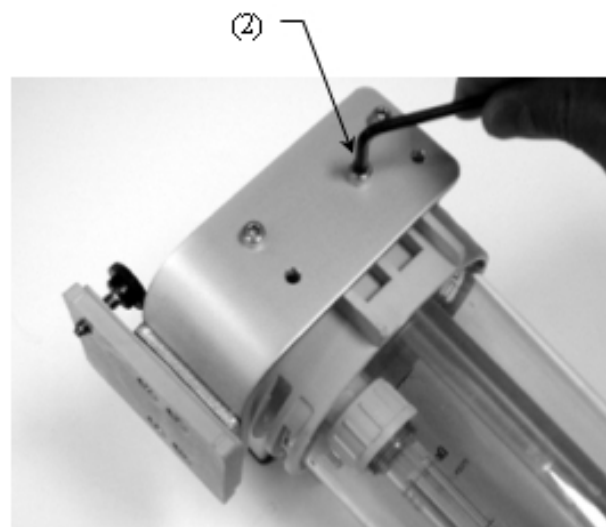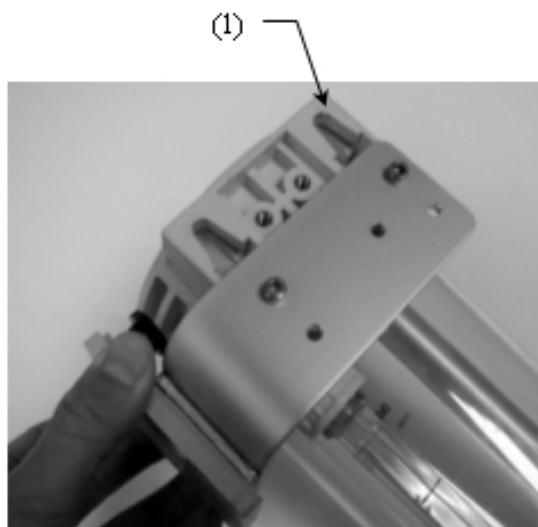

Mounting the Scavenger to Bracket

1. Guide the Slide (Mounting Bracket) into the channel. Move Scavenger to desired height and tighten black Knob.

Fabius GS Configurations



Left Side Attachment



Rear Attachment



Right Side Attachment


FABIUS TIRO ASSEMBLY INSTRUCTIONS:

Configuring the Mounting Bracket for the Application


1. Using a Philips screwdriver, remove the Slide (item 1) by removing 4 screws from the existing configuration and reattach the Slide in the new configuration as shown below.
2. Using the 5/32" hex key (provided), remove the two (2) mounting posts (item 2) from the existing configuration and reattach them in the new configuration as shown below. Ensure that the mounting posts are located in the upper two mounting holes on the bracket.

Mounting Bracket Configurations for Fabius Tiro

Left Side Attachment

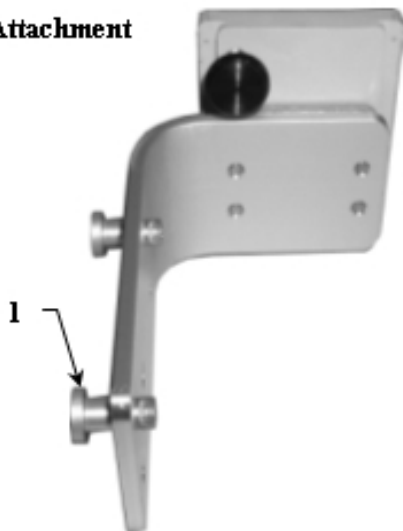



Right Side Attachment

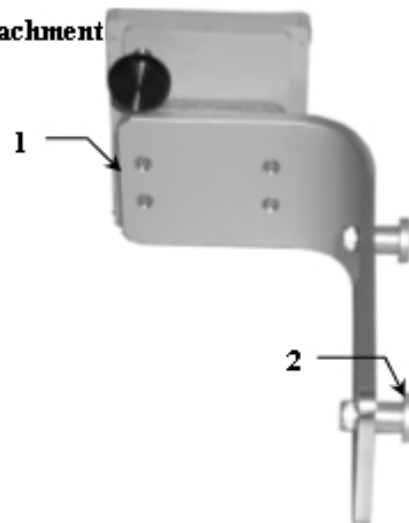

Attaching the Scavenger to the Mounting Bracket for Fabius Tiro

1. Guide the two (2) mounting posts into the mounting slots (item 1) on the rear of the Scavenger.
2. Using the 5/32" hex key (provided), secure Scavenger to Mounting Bracket with the M5 x 18mm socket head cap screw (SHCS) (item 2).

SP00284
9/18/03

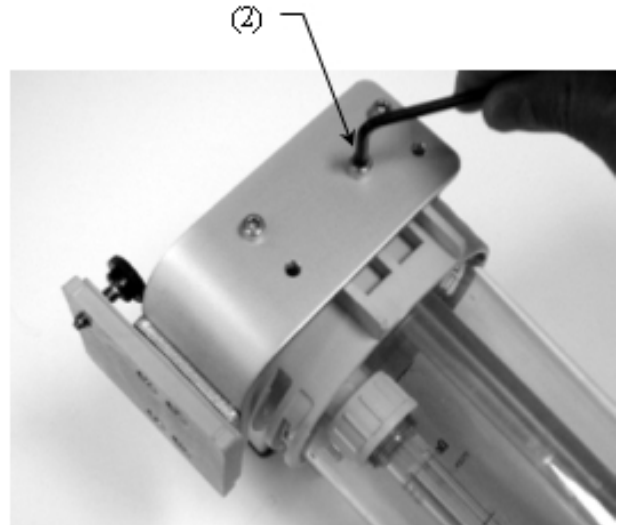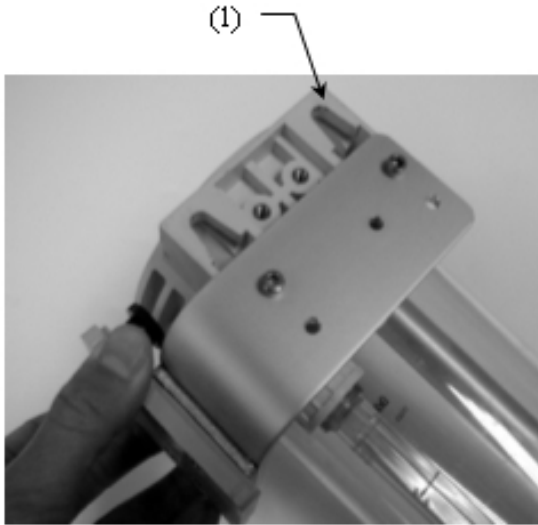

Mounting the Scavenger with Mounting Bracket to Fabius Tiro

1. Guide the Slide (Mounting Bracket) into the channel. Move Scavenger to desired height and tighten black Knob.



Left Side Attachment



Right Side Attachment

Attaching the Scavenger to the Draeger Clamp for Fabius Tiro

1. Align the holes on the Draeger Clamp with the holes in the AGS Scavenger.
2. Thread the two (2) M5 screws with washers through the clamp.



Mounting the Scavenger with Draeger Clamp to Fabius Tiro

1. Guide the Draeger Clamp on the Draeger rail. Move Scavenger to desired location and clamp with black lever.



Draeger Rail Attachment (for use with Right Mount Cosy)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of
Draeger Medical, Inc.
3122 Commerce Drive
Telford, PA 18969
Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5784
Web: www.draegermedical.com
Printed in the U.S.A.